

Using Lab Equipment

Different experiments require different types of equipment. But even though experiments differ, the ways in which the equipment is used are the same.

Beakers

- Use beakers for holding and pouring liquids.
- Do not use a beaker to measure the volume of a liquid. Use a graduated cylinder instead. (See page R16.)
- Use a beaker that holds about twice as much liquid as you need. For example, if you need 100 milliliters of water, you should use a 200- or 250-milliliter beaker.




Test Tubes

- Use test tubes to hold small amounts of substances.
- Do not use a test tube to measure the volume of a liquid.
- Use a test tube when heating a substance over a flame. Aim the mouth of the tube away from yourself and other people.
- Liquids easily spill or splash from test tubes, so it is important to use only small amounts of liquids.

Test-Tube Holder

- Use a test-tube holder when heating a substance in a test tube.
- Use a test-tube holder if the substance in a test tube is dangerous to touch.
- Make sure the test-tube holder tightly grips the test tube so that the test tube will not slide out of the holder.
- Make sure that the test-tube holder is above the surface of the substance in the test tube so that you can observe the substance.

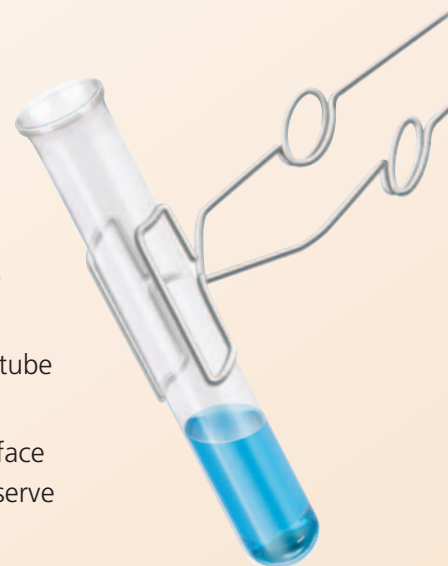

Test-Tube Rack

- Use a test-tube rack to organize test tubes before, during, and after an experiment.
- Use a test-tube rack to keep test tubes upright so that they do not fall over and spill their contents.
- Use a test-tube rack that is the correct size for the test tubes that you are using. If the rack is too small, a test tube may become stuck. If the rack is too large, a test tube may lean over, and some of its contents may spill or splash.




Forceps

- Use forceps when you need to pick up or hold a very small object that should not be touched with your hands.
- Do not use forceps to hold anything over a flame, because forceps are not long enough to keep your hand safely away from the flame. Plastic forceps will melt, and metal forceps will conduct heat and burn your hand.



Hot Plate

- Use a hot plate when a substance needs to be kept warmer than room temperature for a long period of time.
- Use a hot plate instead of a Bunsen burner or a candle when you need to carefully control temperature.
- Do not use a hot plate when a substance needs to be burned in an experiment.
- Always use “hot hands” safety mitts or oven mitts when handling anything that has been heated on a hot plate.

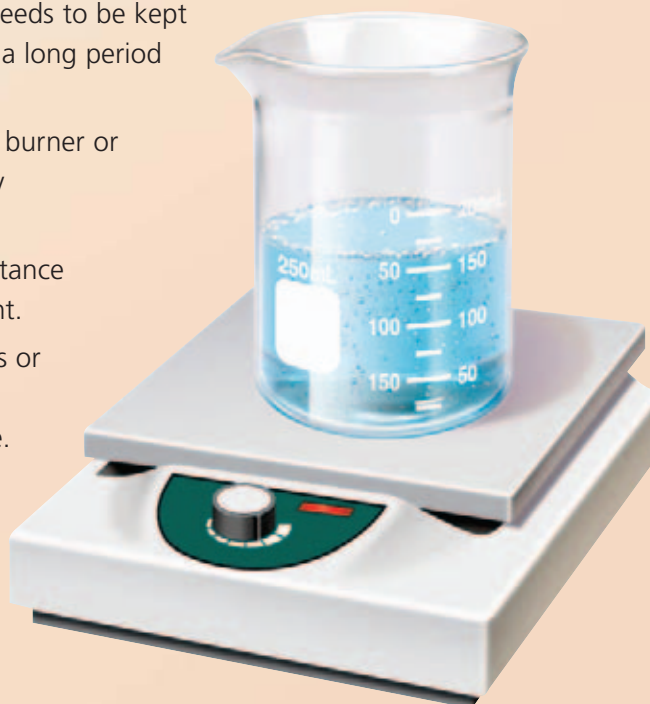

Microscope

Scientists use microscopes to see very small objects that cannot easily be seen with the eye alone. A microscope magnifies the image of an object so that small details may be observed. A microscope that you may use can magnify an object 400 times—the object will appear 400 times larger than its actual size.

Body The body separates the lens in the eyepiece from the objective lenses below.

Nosepiece The nosepiece holds the objective lenses above the stage and rotates so that all lenses may be used.

High-Power Objective Lens This is the largest lens on the nosepiece. It magnifies an image approximately 40 times.

Stage The stage supports the object being viewed.

Diaphragm The diaphragm is used to adjust the amount of light passing through the slide and into an objective lens.

Mirror or Light Source Some microscopes use light that is reflected through the stage by a mirror. Other microscopes have their own light sources.

Eyepiece Objects are viewed through the eyepiece. The eyepiece contains a lens that commonly magnifies an image 10 times.

Coarse Adjustment This knob is used to focus the image of an object when it is viewed through the low-power lens.


Fine Adjustment This knob is used to focus the image of an object when it is viewed through the high-power lens.

Low-Power Objective Lens This is the smallest lens on the nosepiece. It magnifies an image approximately 10 times.

Arm The arm supports the body above the stage. Always carry a microscope by the arm and base.

Stage Clip The stage clip holds a slide in place on the stage.


Base The base supports the microscope.

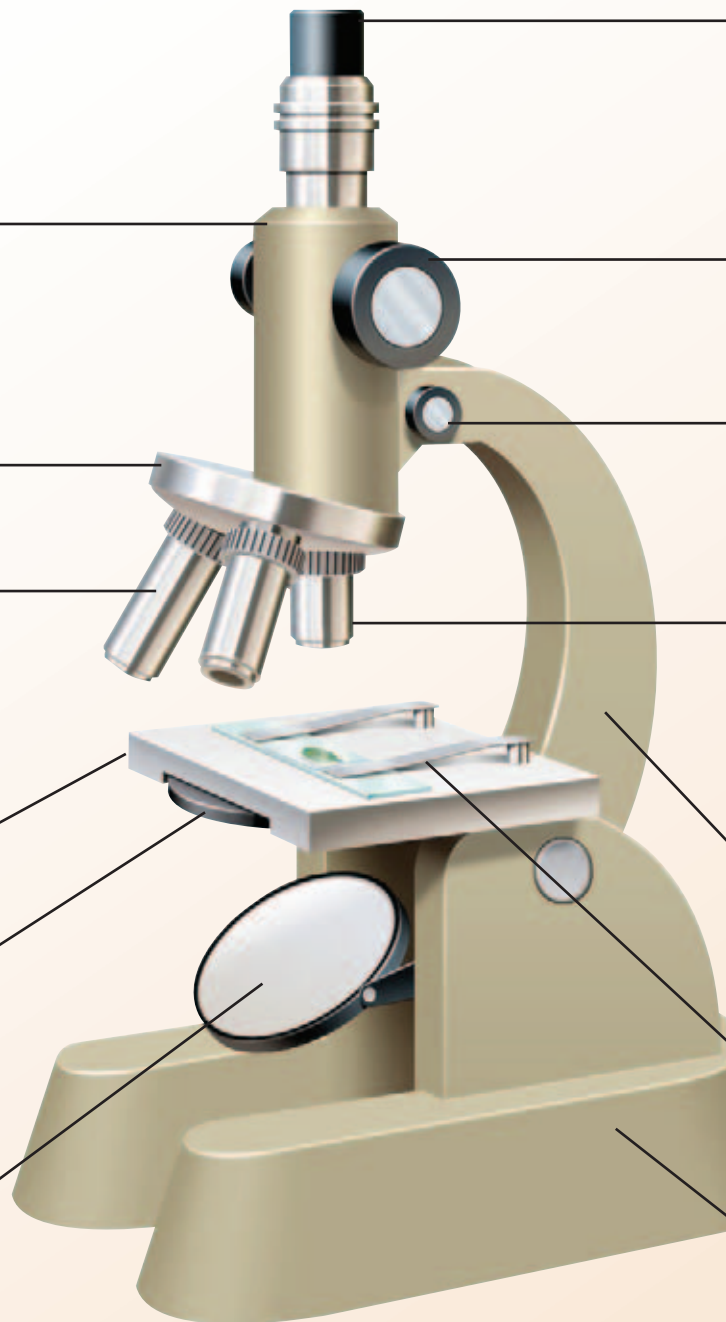

VIEWING AN OBJECT

1. Use the coarse adjustment knob to raise the body tube.
2. Adjust the diaphragm so that you can see a bright circle of light through the eyepiece.
3. Place the object or slide on the stage. Be sure that it is centered over the hole in the stage.
4. Turn the nosepiece to click the low-power lens into place.
5. Using the coarse adjustment knob, slowly lower the lens and focus on the specimen being viewed. Be sure not to touch the slide or object with the lens.
6. When switching from the low-power lens to the high-power lens, first raise the body tube with the coarse adjustment knob so that the high-power lens will not hit the slide.
7. Turn the nosepiece to click the high-power lens into place.
8. Use the fine adjustment knob to focus on the specimen being viewed. Again, be sure not to touch the slide or object with the lens.

MAKING A SLIDE, OR WET MOUNT

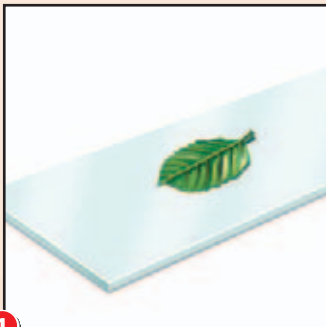

- 1** Place the specimen in the center of a clean slide.



- 2** Place a drop of water on the specimen.



- 3** Place a cover slip on the slide. Put one edge of the cover slip into the drop of water and slowly lower it over the specimen.




- 4** Remove any air bubbles from under the cover slip by gently tapping the cover slip.



- 5** Dry any excess water before placing the slide on the microscope stage for viewing.


Spring Scale (Force Meter)

- Use a spring scale to measure a force pulling on the scale.
- Use a spring scale to measure the force of gravity exerted on an object by Earth.
- To measure a force accurately, a spring scale must be zeroed before it is used. The scale is zeroed when no weight is attached and the indicator is positioned at zero.
- Do not attach a weight that is either too heavy or too light to a spring scale. A weight that is too heavy could break the scale or exert too great a force for the scale to measure. A weight that is too light may not exert enough force to be measured accurately.

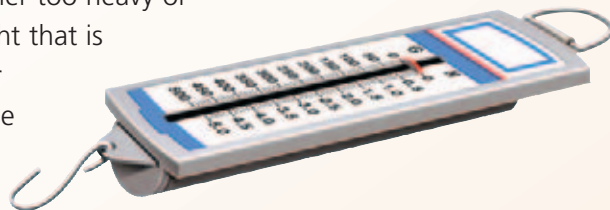

Graduated Cylinder

- Use a graduated cylinder to measure the volume of a liquid.
- Be sure that the graduated cylinder is on a flat surface so that your measurement will be accurate.
- When reading the scale on a graduated cylinder, be sure to have your eyes at the level of the surface of the liquid.
- The surface of the liquid will be curved in the graduated cylinder. Read the volume of the liquid at the bottom of the curve, or meniscus (muh-NIHS-kuhs).
- You can use a graduated cylinder to find the volume of a solid object by measuring the increase in a liquid's level after you add the object to the cylinder.

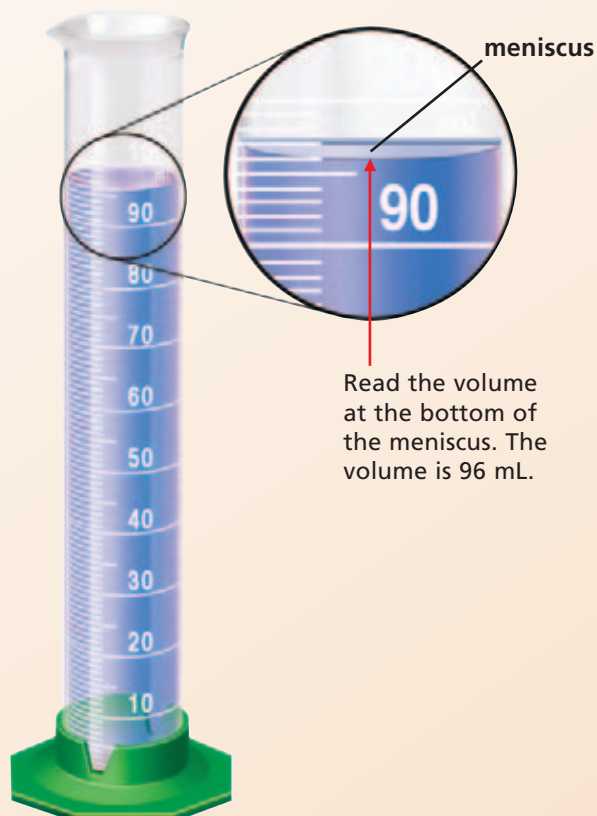

Metric Rulers

- Use metric rulers or meter sticks to measure objects' lengths.
- Do not measure an object from the end of a metric ruler or meter stick, because the end is often imperfect. Instead, measure from the 1-centimeter mark, but remember to subtract a centimeter from the apparent measurement.
- Estimate any lengths that extend between marked units. For example, if a meter stick shows centimeters but not millimeters, you can estimate the length that an object extends between centimeter marks to measure it to the nearest millimeter.
- **Controlling Variables** If you are taking repeated measurements, always measure from the same point each time. For example, if you're measuring how high two different balls bounce when dropped from the same height, measure both bounces at the same point on the balls—either the top or the bottom. Do not measure at the top of one ball and the bottom of the other.

EXAMPLE

How to Measure a Leaf

1. Lay a ruler flat on top of the leaf so that the 1-centimeter mark lines up with one end. Make sure the ruler and the leaf do not move between the time you line them up and the time you take the measurement.
2. Look straight down on the ruler so that you can see exactly how the marks line up with the other end of the leaf.
3. Estimate the length by which the leaf extends beyond a marking. For example, the leaf below extends about halfway between the 4.2-centimeter and 4.3-centimeter marks, so the apparent measurement is about 4.25 centimeters.
4. Remember to subtract 1 centimeter from your apparent measurement, since you started at the 1-centimeter mark on the ruler and not at the end. The leaf is about 3.25 centimeters long ($4.25 \text{ cm} - 1 \text{ cm} = 3.25 \text{ cm}$).




Triple-Beam Balance

This balance has a pan and three beams with sliding masses, called riders. At one end of the beams is a pointer that indicates whether the mass on the pan is equal to the masses shown on the beams.

1. Make sure the balance is zeroed before measuring the mass of an object. The balance is zeroed if the pointer is at zero when nothing is on the pan and the riders are at their zero points. Use the adjustment knob at the base of the balance to zero it.
2. Place the object to be measured on the pan.
3. Move the riders one notch at a time away from the pan. Begin with the largest rider. If moving the largest rider one notch brings the pointer below zero, begin measuring the mass of the object with the next smaller rider.
4. Change the positions of the riders until they balance the mass on the pan and the pointer is at zero. Then add the readings from the three beams to determine the mass of the object.


300 g	position of largest rider
90 g	position of middle rider
+ 3 g	position of smallest rider
<u>393 g</u>	mass of beaker

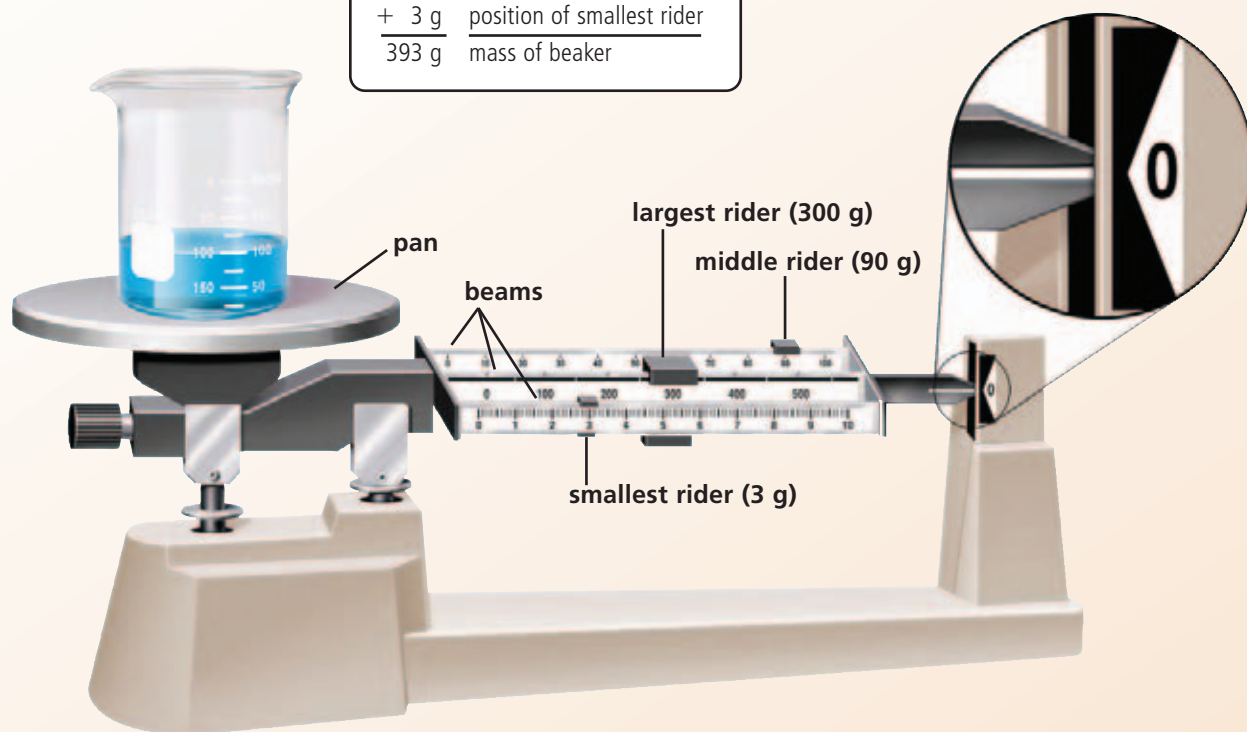

Double-Pan Balance

This type of balance has two pans. Between the pans is a pointer that indicates whether the masses on the pans are equal.

1. Make sure the balance is zeroed before measuring the mass of an object. The balance is zeroed if the pointer is at zero when there is nothing on either of the pans. Many double-pan balances have sliding knobs that can be used to zero them.
2. Place the object to be measured on one of the pans.
3. Begin adding standard masses to the other pan. Begin with the largest standard mass. If this adds too much mass to the balance, begin measuring the mass of the object with the next smaller standard mass.
4. Add standard masses until the masses on both pans are balanced and the pointer is at zero. Then add the standard masses together to determine the mass of the object being measured.

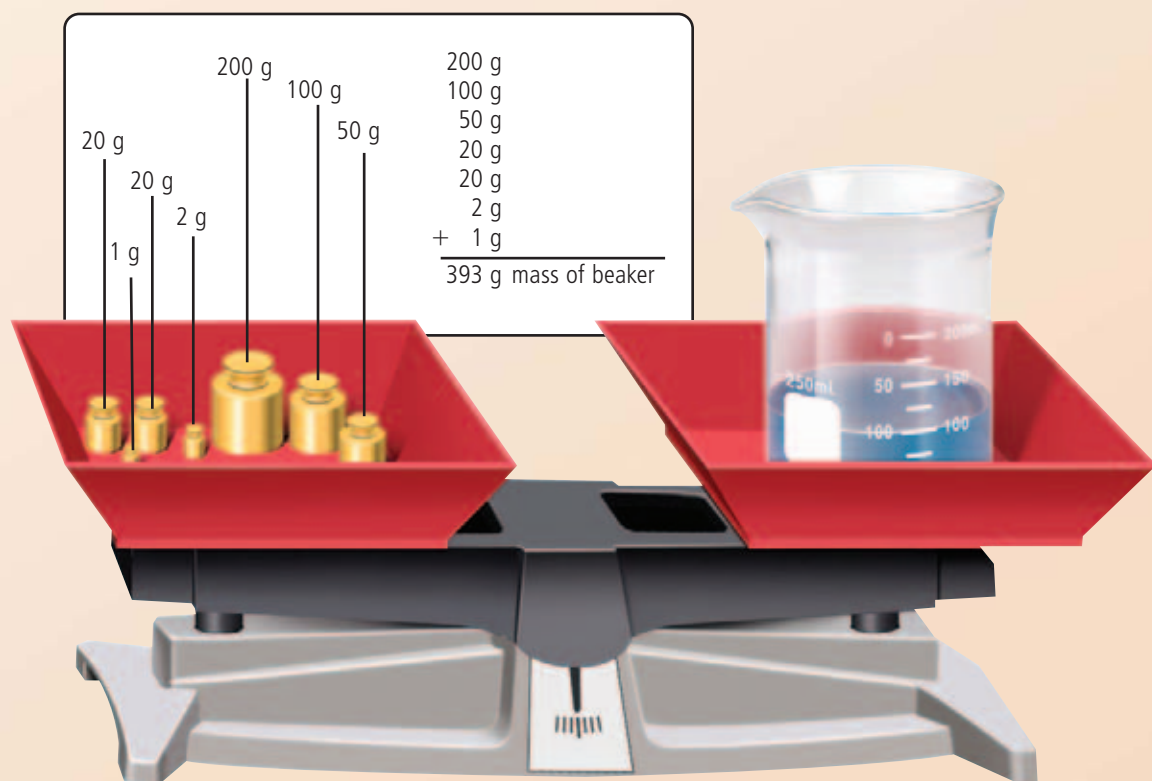

Never place chemicals or liquids directly on a pan. Instead, use the following procedure:

- 1 Determine the mass of an empty container, such as a beaker.
- 2 Pour the substance into the container, and measure the total mass of the substance and the container.
- 3 Subtract the mass of the empty container from the total mass to find the mass of the substance.